


SPLENDOR® UNIVERSAL



PATENT: PCT/BR2018/050285

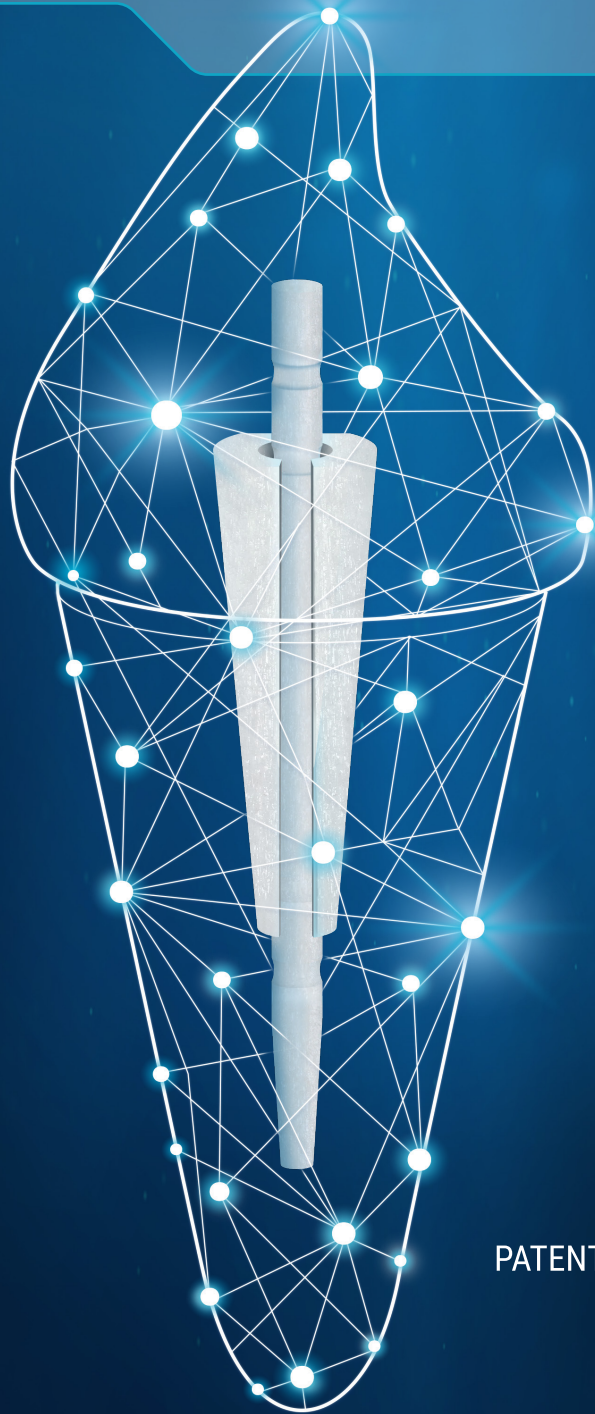

Launch

SPLENDOR[®], beyond universal, intelligent!


A new concept 2nd generation of fiberglass posts

SPLENDOR[®] UNIVERSAL is a new generation of fiberglass posts with a revolutionary concept and easy-to-use design. Unlike conventional posts that require the wear of the tooth structure for adaptation, SPLENDOR[®] UNIVERSAL allows, through patented technology, the single post self-adjustment in all root canal diameters.

**A UNIVERSAL AND SMART POST
FOR ANY SITUATION!**

INDICATION

- Intracanal support for restorations and fixed prostheses.

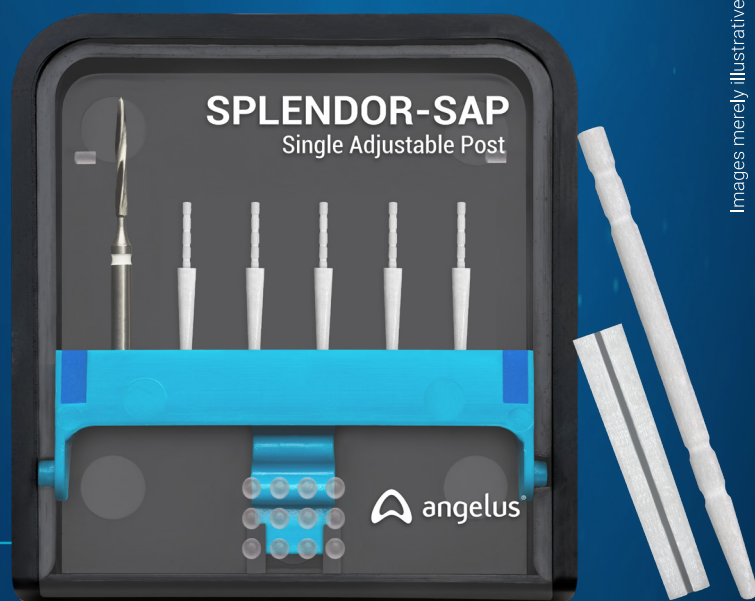

Images merely illustrative.

3 reasons to use SPLENDOR[®]



1. Universal

A single post that adapts to the different root canal diameters (narrow, medium and wide), simplifying the technique and eliminating the need to maintain a stock of posts with different sizes.




2. Self-adjusting

There is no need for drills because it self-adjusts to the canal in its original diameter, adhering to the principles of the conservative dentistry.



3. Anatomical

Promotes high retention in the root canal by the juxtaposition to the root canal walls, with excellent adaptation to the cervical area.

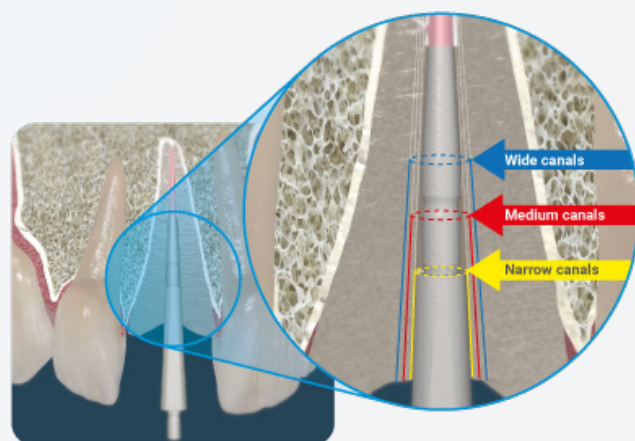

Splendor® in Evidence


Scientific studies performed by various authors demonstrate better outcomes with the use of Universal Splendor:

- Compared to conventional posts, it has a better fit in cervical regions, higher push-out bond strength values, and thus greater retention in the canal.
- Compared to anatomized posts, it provides lower risks of catastrophic tooth root fractures.


For a more comprehensive reading of the articles and further information, please access our interactive content via the QR Code on the side.




Technique of use

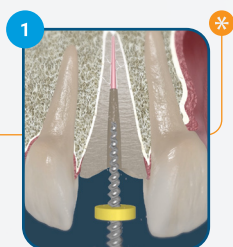

1 Clean the root canal walls with appropriate instruments.

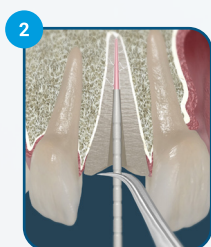


2 Position the post in the canal.

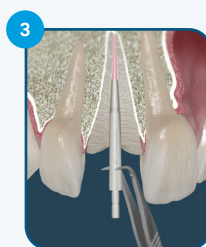


3 Insert the sleeve onto the post and position it as apically as possible.



4 Cut the set - post and sleeve - to the defined height.

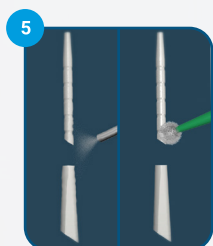

5 Remove the set, clean them separately with alcohol, and dry with air jets. Silanize and leave to dry for 1 minute.

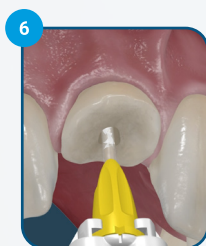

6 Prepare the root canal and insert the cement according to the manufacturer's instructions.



7 First, position the post in the canal and then the sleeve over it.



8 Make the core of the crown restoration.

* The Splendor drill provided in the available presentations is intended exclusively for cleaning the canal.

Case Report

Case provided by **Prof. Leonardo Cunha** and **Prof. Ubiracy Gaião**



1. Testing the SPlendor® UNIVERSAL post in the canal.



2. Positioning the sleeve over the post.



3. Building the core of the crown restoration.

Simplified!

Customer Service

☎ +55 43 9 9919 0244

0800 727 3201

+55 43 2101 3200



/angelusodonto
/angelusodontoglobal
/angelusLATAM



/AngelusOdonto
/AngelusOdontoGlobal
/AngelusOdontoLATAM

www.angelus.ind.br

